



The next generation of  
neuromuscular diagnostics





## Ultra-Simple. Ultra-Certain.

When it comes to EMG, NCS and EP testing, the best tool is one you know you can rely on. Building on the success of the VikingQuest and previous UltraPro models, Natus is proud to introduce the next generation UltraPro S100. Developed in close collaboration with physicians and technologists, the UltraPro's easy-to-use, familiar interface and superior signal quality help assess patients with a greater level of certainty.


Enhance performance with the UltraPro:

- Simple setup to start performing studies faster
- Familiar, customizable workflows and user-friendly interface
- Superior signal quality for High-Definition EMG™ performance
- Automatic selection of AANEM reference values
- Added and improved testing capabilities

HD EMG™



Includes electrical stimulator probe

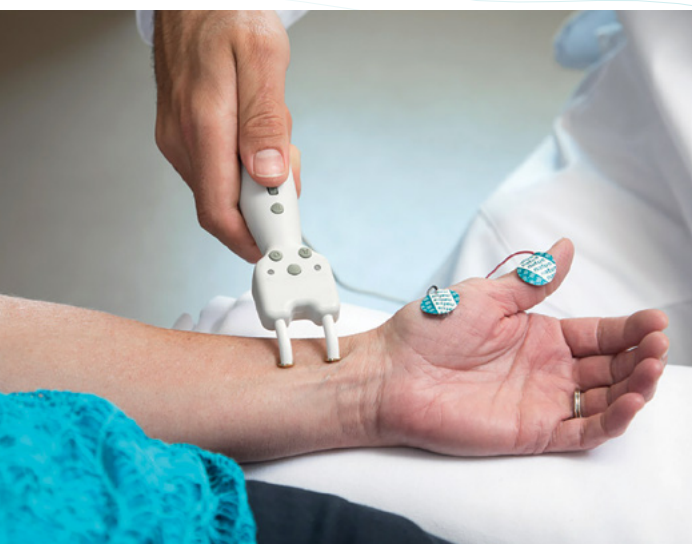

## The familiar way to reach an informed diagnosis

### Faster, simplified workflow from patient to report

Designed with the user in mind, the UltraPro's enhancements in hardware and software make it easier for physicians and technologists to get better data when testing for muscle and nerve disorders — and give patients better answers.

4- channel amplifier option available



### A system improved. Easy to use.

- Built-in control panel with integrated speakers; electrical, visual and auditory capabilities; choice of 3- or 4-channel amplifier; connectors for optional devices; and laptop computer with system software
- Easy-to-use interface with intuitive control panel buttons and color-coded function keys — already familiar to Natus EMG legacy users
- Programmable stimulator buttons for key function control at your fingertips
- Triple-programmable footswitch for hands-free operation
- Simple, three-cable setup to begin testing sooner
- Easily connects to peripheral devices

### Trusted software. Reliable data.

- Built with industry-leading Natus Elite™ software powered by Viking or Synergy, with full patient and protocol compatibility with existing Natus products
- Simplified, customizable workflows speed up frequent operations like rollback, roll forward and edit average to generate reports faster
- Includes AANEM reference values, automatically selected based on patient age, height, gender and body mass index
- High-Definition EMG™ provides a resolution of 4,800 points per trace for saved data in routine motor and sensory nerve conduction studies produces high-quality signals
- Natus Elite software with HD EMG™ includes additional, automated testing capabilities and improved integration to keep patient information easily accessible — and more secure



User-friendly, color-coded controls

### Comfortable for providers and patients.

- Available in two convenient form factors to accommodate your needs: a lightweight system that travels easily in roller case or mobile, cart-based system
- Upgraded control panel allows care providers — and controls — to remain right at the patients' side
- Faster workflows and minimal repeat stimulation reduce time and discomfort of testing procedures
- Add confidence to your EMG testing with InVisus Neuromuscular Ultrasound\*

\*For InVisus availability please contact your local representative



## Advancing neuromuscular diagnostics. Together.

By choosing the UltraPro, you're choosing more than an enhanced system with better signal quality — you're choosing a partnership with leaders in EMG diagnostic development, backed by over 60 years of clinical knowledge. Natus is committed to supporting our customers and their teams to ensure a seamless transition. Customers can rely on Natus for more than superior technology.

- Access to our full-time clinical specialist support team to provide educational sessions
- Large field service team and distribution partners for installation assistance, rapid, on-site troubleshooting and technical support
- A robust help menu and extensive video training library for assistance with hardware, workflows, diagnostic tests, system configurations and more; also available off-site via the Neuro Training Academy. Membership is free and open to all. Join today at [neuro-training.academy](http://neuro-training.academy)



## Natus Support

At Natus, we strive for excellence in customer and technical service.

### Here's how we can help:

- Accessible and effective technical support
- Definitive technical documentation and knowledgeable installation teams
- Optional extended warranty and service coverage programs
- Comprehensive, flexible customer training courses

## Supplies

Convenient, complete, trusted.

Natus supports the full spectrum of EMG care, providing a complete portfolio of EMG supplies for a seamless solution.

- Dedicated and knowledgeable customer support
- Streamlined order processing
- Convenient online ordering (US Customers only)  
[Natus Medical Store – natusmedicalstore.com](http://Natus Medical Store – natusmedicalstore.com)

To learn more about Natus products, contact your local distributor or sales representative.

US Customers Call: **1-800-356-0007**

International Customers Call: **+1-608-829-8500**

## Healthcare solutions with one thing in mind. You.

©2021 Natus Medical Incorporated. All Rights Reserved. All product names appearing on this document are trademarks or registered trademarks owned, licensed to, promoted or distributed by Natus Medical Incorporated, its subsidiaries or affiliates. 019824 RevD1

# natus®

Natus Medical Incorporated

[natus.com/ultrapro](http://natus.com/ultrapro)